

New

Vial feeding pox model

Vial feeding fullpox™ 185 syringe

The fullpox™ 185 syringe is a robust, dual-needle tool specifically designed for effective delivery of fowl pox vaccination. It is built for durability and ease of use, meeting the needs of veterinary professionals in poultry farming.

- Natural ergonomics, effortless activation
- Retractable dual-needle solution tailored for fowl pox vaccination
- Economical, fully interchangeable parts for cost-effective maintenance and long-term use
- Matching expected quality level from Socorex
- Fully autoclavable 121°C /250°F
- Re-usable packaging keeping syringe protected



Comfortable by design

Ergonomic handle and smooth activation grant a substantial working comfort.



Dual needle

The syringe delivers precise 0,01 mL doses with every inoculation, ensuring accuracy for effective fowl pox protection. Its dual eye-of-the-needle design supports efficient and reliable intradermal delivery, while retractable needles enhance operator safety by

minimizing the risk of accidental injuries.



Bottle with cap

Specially designed 15 mL bottle ensures up to 1'500 seamless vaccine deliveries. Reusable and made from high-quality polypropylene (PP), both bottle and cap are fully autoclavable.


Maintenance and sustainability

Stainless steel body and handle made from glass fibre-reinforced resin offer excellent robustness. The instrument can be fully disassembled and reassembled making cleaning and maintenance a breeze.

Ordering information

Standard supply includes syringe, bottle, two spare needles, Allen key and quick guide.

| Description | Packaging | Cat. No. |
|--------------------------|-----------|------------|
| Syringe fullpox™ 185 | 1 / pk | 185.1.0010 |
| Spare inoculation needle | 2 / pk | 185.902 |
| Spare inoculation needle | 10 / pk | 185.910 |
| Bottle with cap, 15 mL | 1 / pk | 185.528 |
| Barrel cap | 1 / pk | 185.721 |
| Screw for needle | 1 / pk | 185.801 |

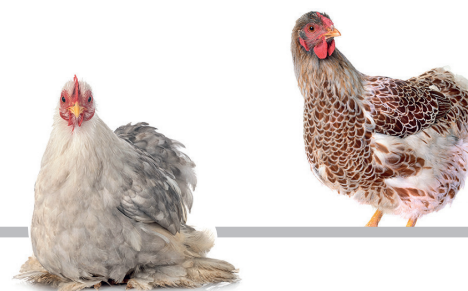

Pro Tip

To maximize vaccination efficiency, pre-fill several bottles with vaccine solution and store in cool conditions for quick access during campaigns.